

Manafam



Blunt Micro-Cannula

Product Catalogue

www.tebyanmed.com




Differences *Are Visible* *Tangible* *Effective*




Optimised Tip Shape

Using the UltimateFill technology, the tip of the cannula is well filled and obtains the appropriate shape to move through the layers of the skin smoothly without any damage to the vessels or tissues.

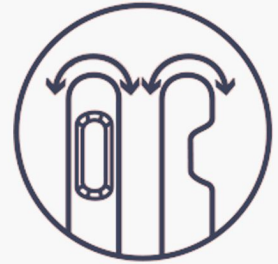


Soft edges in the slot

In creating the slot in this microcannula, MicroPolish technology is used to create smooth and rounded edges on the sides of the slot and eliminate the possibility of any internal bleeding.

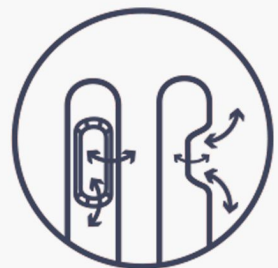


Strong and slippery silicon coating on steel

High-quality silicone, Corning formulation and advanced coating technology have resulted in a very clear and thin coating on the surface of the cannula, making it easy to move while working.

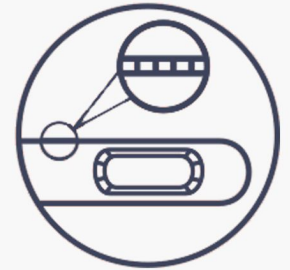


Thinner tube wall and larger internal diameter

The thinner the tube wall, the larger the inner diameter of the cannula. This extra space helps the material move easily through the cannula.

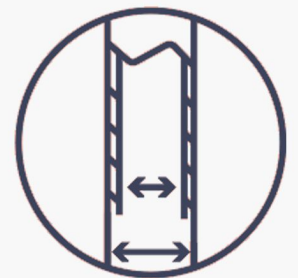


Standardized hub slope, fix and reliable

The plastic base part (hub) of the product has full compliance with the ISO80369 standard to be easily attached to standard syringes and prevent any leakage.

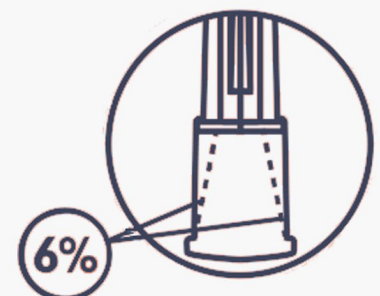

Combining of stiffness and flexibility for perfect performance

The combination of strength and flexibility will help you both in the case of deep injections where you need strength and under the skin where you need flexibility.

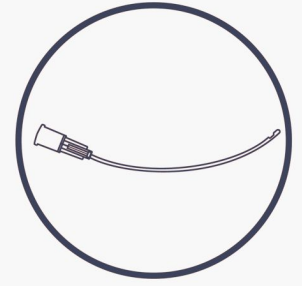

Manafam

First package



1 Microcannula and 1 hypodermic needle of appropriate size

Second package














50 peaces of first package in a box

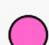Third package



48 box in a carton

Specifications

Size	Color	Diameter (mm)	length (mm)
G16		1.6	90
G16		1.6	110
G18		1.2	70
G18		1.2	90
G18		1.2	110
G21		0.8	70
G22		0.7	50
G22		0.7	70
G23		0.6	50
G25		0.5	50
G27		0.4	50

* Each microcannula is supplied with a hypodermic needle.
The size of the hypodermic needle is related to the size of the microcannula.

Blunt Micro-cannula

Manafam

Produced and registered by



Custom Medical Needles Production
and Metal Tubes Forming

Address: No. 2, North Milad Blvd, Golgoon
Industrial City, Tehran, Iran

Phone: (+98) 21 77 03 86 27

Messangers: (+98) 920 120 2608

www.tebyanmed.com